



URSATEC  
WHEN PRESERVATIVE FREE MATTERS®

SoftBreezer®  
gentle breathe



Ursatec dosage systems  
for preservative free inhalation

# Preservative free for a healthy life®

## The philosophy behind Ursatec Verpackung GmbH


Ursatec Verpackung GmbH is dedicated to support the global health care industry in avoiding preservatives in any kind of product. Health care administrations recommend the development of preservative-free products worldwide.

*„Devices and manufacturing processes shall be designed in such a way as to eliminate or to reduce as far as possible the risk of infection to patients, users and, where applicable, other persons.“*

Interinstitutional File: 2012/0266 (COD); Brussels, 27 June 2016 (OR. en); 10617/16  
Council of the European Union

## SoftBreezer® – innovative spray nozzle for gentle breathe

The SoftBreezer® spray nozzle unit developed by Medspray is based on highly precise orifices in a small surface.

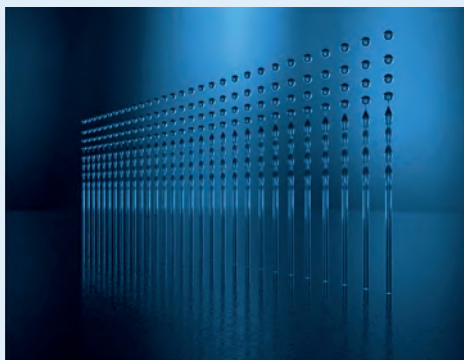

Pressing water through the SoftBreezer® spray nozzle unit results in multiple jets which will automatically break up in small droplets (Rayleigh principle). The resulting spray mist is mono-disperse which means that all droplets have the same size for good lung deposition.

# SoftBreezer® – refresh your life

## Refresh your life

First preservative-free mobile nebulizer for consumer health products such as isotonic saline inhalation solutions. SoftBreezer® is a manually actuated, pocket size nebulizer device for aqueous solutions, releasing a mono-disperse aerosol, targeting the airways very efficiently.

The SoftBreezer® preservative-free spray pump is based on Ursatec's proven spray devices used for over 20 years.

## Key Features

Suitable for all kind of preventive products.

- Gentle breathe
- Preservative free
- Propellant free
- Pocket size
- Proven 3K® technology
- Microbiologically safe

## Data sheet

---

**Pump delivery PD (water placebo)**

25 mg

---

**Airflow at 2 kPa**

10 -20 lpm

---

**Container fill**

5 ml or 10 ml

---

**Number of actuations/container**


200 or 400

---

**Microbial integrity of container in use**

min. 90 days

---

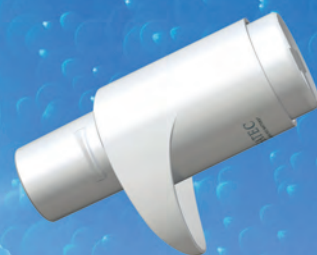

Ursatec Verpackung GmbH  
Schillerstraße 4  
D-66606 St. Wendel

Telefon: 0049 (0) 6851 - 80 26 0  
Telefax: 0049 (0) 6851 - 80 26 14

E-Mail: [info@ursatec.com](mailto:info@ursatec.com)  
Internet: [www.ursatec.com](http://www.ursatec.com)